

# Using machine learning to perform image quality control: A computer exercise using CellProfiler & CellProfiler Analyst software

Beth Cimini, Broad Institute

## Background information:

High-quality image analysis relies on having high-quality images to analyze; the so-called ‘Garbage in, garbage out’ principle. If your image set is small it is straightforward to curate it and remove images that are not in focus or that contain debris that might disrupt segmentation. In a high content screen containing tens or hundreds of thousands of images, this is essentially impossible.

The images in this experiment come from the [Broad Bioimage Benchmark Collection](#). They are 240 of 69,120 fields of U2OS cells treated with a panel of 1600 known bioactive compounds and imaged in five channels for a so-called CellPainting assay- see Gustafsdottir et al, 2013 for more information. Given the large size of this experiment, it was necessary to come up with an automated quality control method.

To learn more about machine learning for image quality control, please see Bray et al 2012 and references therein.

## Goals of this exercise:

This exercise will teach you how to measure aspects of image quality with CellProfiler and how to use the machine learning tools in CPA to create robust rules to distinguish low-quality images. It additionally will show you how to incorporate these rules into an image analysis pipeline so that they may guide your downstream data analysis.


## Materials necessary for this exercise:

These 1200 images (240 sites in 5 channels) represent 120 wells from a single 384 well plate, either mock treated with DMSO or treated with a variety of bioactive compounds.

It is additionally expected that you are generally familiar with the tools in CellProfiler and CellProfiler Analyst, preferably after completing the Translocation tutorial or a similar introductory exercise.

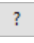
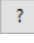
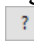
### 1. Start the provided QC pipeline on the BBBC022 dataset

In order to do quality control, we need to first measure the images in many ways. This will allow us to do machine learning to use the measurements to identify the good images from the bad.

- Start CellProfiler by double-clicking the desktop icon 
- Drag and drop the ‘BBBC022\_QC.cppipe’ file into the ‘Analysis modules’ box. Three modules should appear- two copies of ‘MeasureImageQuality’ and one ‘ExportToDatabase’ module.
- Drag and drop the ‘BBBC022\_20585\_AE’ folder into the ‘File list’ box. It should automatically populate.
- Hit the Analyze Images button.
- You will get an error popup – ‘No prior instances of objectgroup were defined. Are you sure you want to continue?’ Hit ‘Yes’.
  - This is due to the fact that there are no objects defined in the pipeline. It will not affect the final results.

- The pipeline should take about 10-15 minutes to run. At the end, three files should be created- a .db database file, a .properties text file, and a .workspace text file.


## 2. Examine the QC pipeline (~15 minutes)

- While the pipeline is running, take some time to look over the pipeline and make sure you understand the various parts. You will probably find it easiest to open a second CellProfiler window and open the pipeline in it.
  - *Metadata*
    - Plate, well, site, and channel information is being extracted from the file and folder names
  - *NamesAndTypes*
    - The images are matched to their channels by metadata and grouped by Plate, Well, and Site
  - *Groups*
    - This module is not used
  - *MeasureImageQuality*
    - ‘Include the image rescaling value?’, ‘Calculate intensity metrics?’, and ‘Calculate saturation metrics?’ are all set to yes. Open the help  button to read about the measurements exported by each of these options.
    - ‘Calculate blur metrics?’ is set to yes and given multiple spatial scales- 2, 5, 10, and 20 pixels. Open the help  button to read about the different methods used to determine whether an image is in focus.
  - *MeasureImageQuality*
    - In this second image ‘MeasureImageQuality’ module, automatic thresholds are calculated for the ‘OrigHoechst’ and ‘OrigPh\_golgi’ images using specific threshold metrics. Read about why by opening the help  button.
      - The choice of these two images and these two thresholding methods was not random; it corresponds to the methods used to segment nuclei and cell bodies in the corresponding analysis pipeline.
  - *ExportToDatabase*
    - This module has 4 main parts, separated by breaks
 
      - **Creation of the database:** these options determine what kind of database will be created (SQLite in this case), whether CellProfiler should overwrite or add to previously existing data, and what the names of the tables will be.
      - **Creation of the CPA properties file:** these options pre-set some parameters for CellProfiler Analyst to help make it easier for you to use such as
        - Setting the plate type to 384, and stating that the Image\_Metadata\_Plate and Image\_Metadata\_Well determine which plate and well to assign each image to
        - What color to assign each channel in the image display tools in CellProfiler Analyst
        - That you are classifying ‘Image’s, not objects
      - **Creation of the CPA workspace file:** instead of having to explore all the parameters at random to find features that separate good images from bad images, these options allow you to auto-create graphs for the features that are likely to be the most important. Here we’ve pre-populated 12 graphs for you-10 scatter plots and 2 histograms.
      - **CellProfiler specific adjustments:** besides determining the location of the output database, most of these options deal with how CellProfiler should treat the objects it identifies (should it calculate means/medians/standard deviations for all of them, should it combine their measurements into a single table, etc). As we are not identifying any objects, we don’t need to worry about these.

## 3. Open the CellProfiler Analyst workspace and determine reasonable parameter cutoffs (~20 minutes)

In the first step of the quality control pipeline, we'll look at graphs of how various measurements are distributed in the population. This allows us to get a sense of what the "normal" values are for each parameter. We'll set gates around the "abnormal" values for later usage. You can feel free to graph and examine additional parameters as well if you are interested in exploring the dataset.

First-pass identification of the abnormal images is important since they are hopefully going to be rare in our overall dataset. By identifying as many abnormal images as possible in this and the next step, we'll hopefully have a large and robust enough training set that we can train a good classifier in later steps.

- Once your pipeline has finished, close CellProfiler and start CellProfiler Analyst by double-clicking the icon  .
- Load the properties file you've just created.
- Go to File-> Load Workspace. 12 graphs should be created – 5 scatter plots about focus (ImageNumber vs ImageQuality\_PowerLogLogSlope), 5 scatter plots about oversaturation (ImageQuality\_StdIntensity vs ImageQuality\_PercentMaximal), and 2 histograms showing the automatic thresholds calculated for the nuclei and cell channels.
- For some or all of the graphs (preferably at least two of each type), examine the shape of the data. Are there obvious outliers present?
  - In the 'gate' dropdown, select CREATE NEW GATE and give it a name (like 'HighSat\_ER' for high outliers in the oversaturation graph of the ER channel).
  - Draw your gate, then right-click (PC)/Command-click (Mac) to bring up a window. You should see two options named 'Show images in gate' (one or both of these will be renamed in the next CPA release for the sake of clarity)
    - The top one will open each of these images in a separate window. This is generally easiest if only a small number of images (<5) are in the gate.
    - The bottom one will open a thumbnail viewer with thumbnails of each image; you can then double click on the thumbnail to open the full image. This is much faster when there are a lot of images inside the gate.

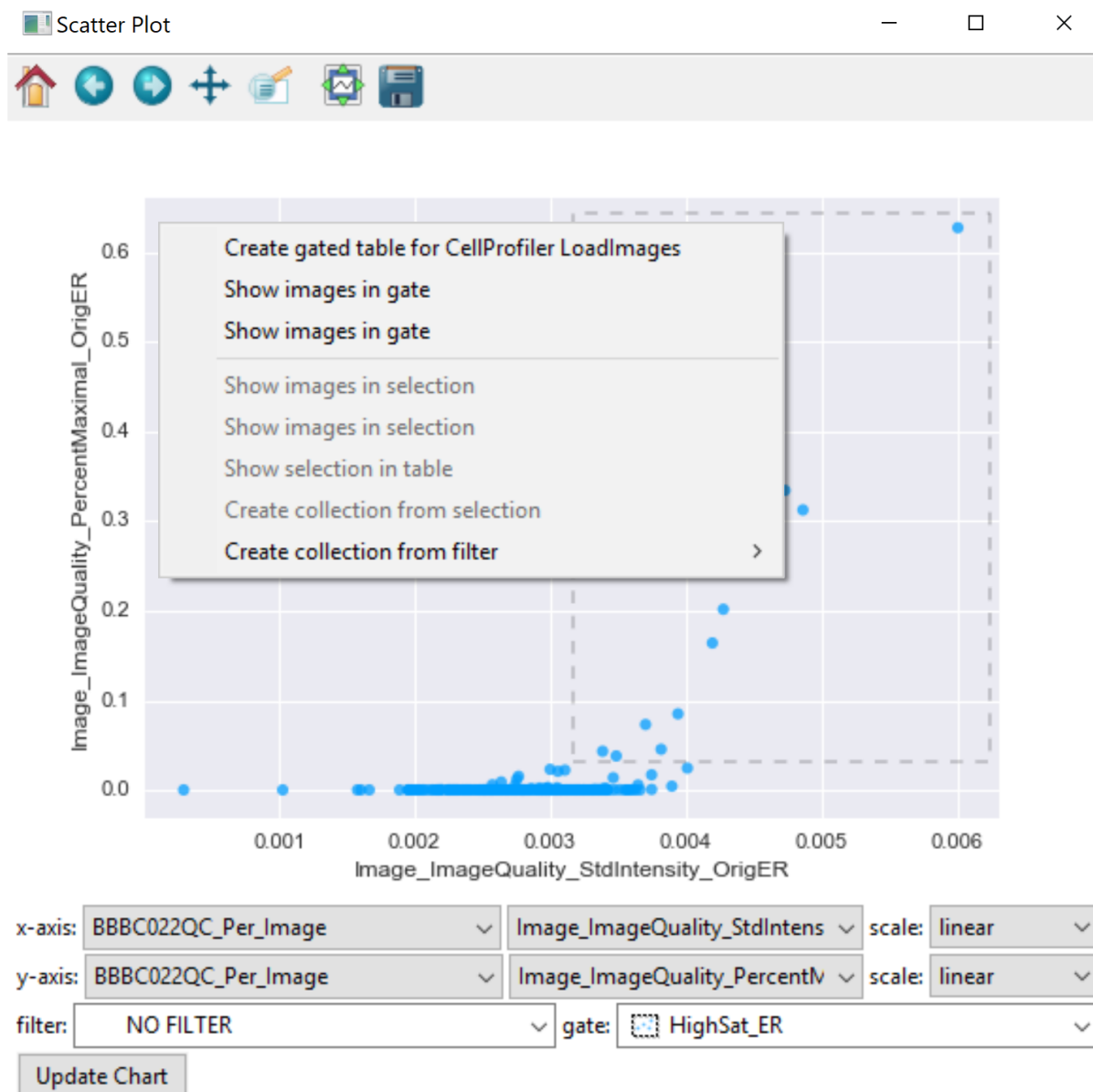


Figure 1: Examining the objects inside a gate.

- Look at some of the images inside the gate ( you may need to click ‘Show controls’ to adjust the zoom, contrast stretch, etc). Are they in focus, and are they free of bright debris? If they seem all or mostly like they’re high quality, adjust the gate to look at a different/smaller subset or delete it altogether by using the Gate Inspector (‘gate’-> ‘MANAGE GATES’).

#### 4. Optional- use the PlateViewer tool to check for other features to gate on (~10 minutes)

If you want to see if you can find additional features that might distinguish good images from bad images, feel free to explore the feature set more thoroughly.

- The PlateViewer tool and either choose a few ImageQuality measurements at random or just use the arrow keys to scroll up and down through the measurements.
  - A good candidate measurement will be relatively constant within most of the wells but have a small number of outliers (Fig 2).
  - Once you have a good candidate, double click on an outlier well to open the images from that well. If one or both seem very aberrant, you may wish to look at this measurement further.

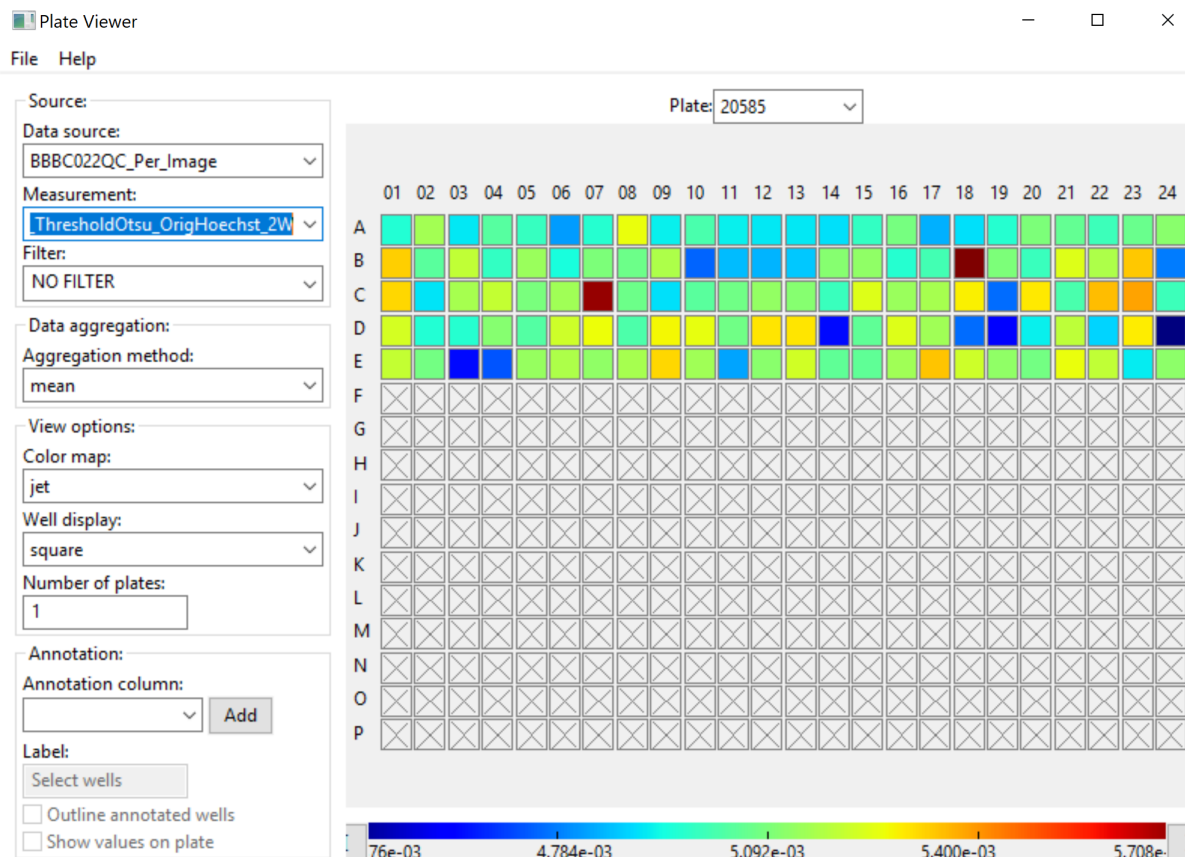


Figure 2: Using the PlateViewer tool to find features that may be good markers of poor quality images.

- Once you have a measurement you'd like to gate based on, open that Histogram or the Scatter Plot tools and graph all the cells based on that (if you're using the Scatter Plot tool, you may wish to graph your chosen feature on the Y-axis vs ImageNumber on the X-axis as in the Workspace plots, but feel free to play around). If you find a gate that seems logical to make, proceed as in Step 3.

## 5. Create filters based on the cutoffs you've determined (~10 minutes)

Now that we've created gates around our poor quality images, we need to convert them into filters so that we can access them in the Classifier tool. This process will be unnecessary in later versions of CellProfiler Analyst, but knowing how to manipulate the properties file is a worthwhile skill.

- When you feel like you've examined enough graphs, it's time to create filters based on the gates you've determined contain bad images. There are two ways to do this:

- **Edit the properties file manually (recommended):**

This method has the advantage that the filters you create are permanent; they won't be lost when you close and re-open CPA. It does however require you to edit the properties file, and potentially requires you to troubleshoot what's wrong with the SQL call if you type it in incorrectly.

- From whatever graph you happen to have open, open the Gate Inspector by calling 'gate'-> 'MANAGE GATES'.
- Open the properties file in a text editor.
- In the Image Filters section, create a new line and begin it with

```
filter_SQL_Blurry = SELECT ImageNumber FROM BBBC022QC_Per_Image WHERE
```

- For each of the blur or threshold gates in your Gate Inspector, copy the text followed by the word OR; as an example if your first gate was:

```
BBBC022QC_Per_Image.Image_ImageQuality_ThresholdOtsu_OrigHoechst_2W BETWEEN
"0.0020927971719" AND "0.0042271321255"
```

Edit your SQL statement to read

```
filter_SQL_Blurry = SELECT ImageNumber FROM BBBC022QC_Per_Image WHERE
BBBC022QC_Per_Image.Image_ImageQuality_ThresholdOtsu_OrigHoechst_2W BETWEEN
"0.0020927971719" AND "0.0042271321255" OR
```

DO NOT hit enter or otherwise add new lines between the different measurements; the statement must be one continuous line in order to work.

- Omit the 'OR' after your last gate measurement.
- Repeat for the saturated gates, again omitting the last OR, following the statement

```
filter_SQL_Saturated = SELECT ImageNumber FROM BBBC022QC_Per_Image WHERE
```

- Save and close the properties file, then re-open CPA.

#### ○ Create the filters inside CPA (not recommended)

This method does not require you to manually edit the properties file, but when filtering based on multiple features it takes longer to set up, cannot be edited, and would need to be repeated every time you closed and re-opened CellProfiler Analyst. It is therefore not recommended for use in this tutorial, but is presented in case you cannot or will not edit the properties file.

- From whatever graph you happen to have open, open the Gate Inspector by calling 'gate' -> 'MANAGE GATES'.
- Copy the text of the relevant gates into a text editor so that you can access it later.
- Determine which is the important cutoff for each gate (ie for a gate where values are too high, it's the minimum value, for gates on values that are too low it's the maximum value)
- Open the filter creation dialog by selecting 'filter' -> 'CREATE NEW FILTER'
- Name your filter Blurry
- Enter each important cutoff from the blur or threshold gates, selecting 'Add column' after each and changing the separator statement to OR
- Create a filter named Saturated containing the important cutoffs from the saturation gates

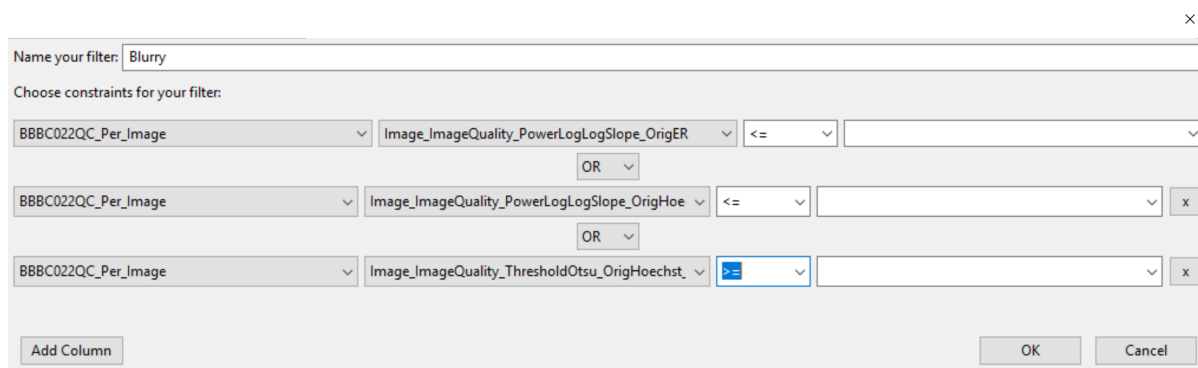


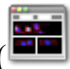
Figure 3: Creating filters inside CPA

## 6. Create classifier rules to distinguish good from bad images (~30 minutes)




Creating gates based on the 12 measurements we graphed has helped us identify some low quality images so far, but we are not utilizing the rest of our measurements and we can't be certain we've found them all. By creating a machine learning classifier we can find features that correlate well with bad images in our experiment and increase our likelihood of identifying additional low quality images.



- Open the Classifier tool ( , 2<sup>nd</sup> on left).
- Fetch 50 random images from 'experiment', then 25 each from 'Blurry' and 'Saturated'
- Create a good classifier by dragging images to the positive and negative bins and training with the Fast Gentle Boosting classifier. If needed feel free to refer to the Translocation tutorial instructions for reminders of how to fetch more images, check the confusion matrix of your classifier, etc.
  - In a larger experiment, you might decide to add more bins (for images that are out of focus, have debris, contain the well edge, etc), but stick to two bins for this experiment.
  - Don't only fetch additional images from your filtered sets, make sure to pull some from 'Random' as well, since your gates may not contain every bad image in the experiment and since you want a robust set of good-quality images for your positive bin.
  - You will likely have to use more than 5 rules to get good classification- try setting this number to 25
- When you're happy with your classifier (this is a small experiment, so perfect classification will not be possible, but aim for >60% accuracy in each bin at least), examine the rules created by using 'Advanced' -> 'Edit rules...'
- Remove any rules that DO NOT have to do with ImageQuality (like 'GroupIndex, etc), then copy the remainder into a text editor and save as a '.txt' file. Keep the text file open for now.
- Save your training set for future reference if desired, then close CellProfiler Analyst.

## 7. Add quality control steps to an analysis pipeline(~15 minutes)

If you have time, you can add the list of rules you identified in your machine learning classifier to the CellProfiler pipeline that corresponds to this data set. This can be valuable to do in processing your own large datasets, especially if time or computational power is at a premium.

- Reopen CellProfiler, then open the BBBC022\_Analysis pipeline (provided or available as part of the "AdvancedSegmentation" tutorial).
- Reopen the BBBC022\_QC pipeline for reference, either in a separate CellProfiler window or in a text editor.
- Click on the  button to open the 'Add modules' window
- If your rule set contains any rules based on 'ImageQuality\_Threshold', add two 'MeasureImageQuality' modules (from the 'Measurement' category) to the very top of your pipeline, otherwise add one.
- Configure the 'MeasureImageQuality' module(s) so that they'll provide all the measurements specified in your rules file; you can refer to the settings in the QC pipeline and/or the rules file to determine which these are.
  - For measurements that end in a number, like 'Image\_ImageQuality\_Correlation\_OrigER\_10', the number refers to the 'Spatial scale for blur measurements' scale specified. You started by measuring the blur at 4 scales, but you may or may not need all of them based on the rules selected.
- Add the 'FlagImage' module from the 'Data Tools' category.
- Leave the flag category as 'Metadata', and the name as 'QCFlag'
  - If you had created many bins in your classification, you would set this name as something more descriptive, such as 'QCDebris', since you'd have a different flag name for each bin in your classifier.
- Leave the 'Skip image set if flagged?' setting at no
  - Your classifier probably wasn't incredibly accurate, so setting this to 'yes' might result in throwing out real data while still not excluding all the bad images. Nevertheless, if you had far more images than you felt you needed you *could* set this to yes.
- Set 'Flag is based on' to 'Rules' and provide the path and file information to your rules text file.

- Set the Class number to whatever class in your rule set represents the bad images- this is 2 if you kept the default bins and dragged bad images to the 'Negative' bin.
- Save your pipeline.